



Yes, you can.®

rea  **Azalea**® /  **Azalea**® Assist

**Offers excellent seating comfort and many adjustment possibilities!**

The Rea™ **Azalea**® is the natural choice for the moderately active user as well as for the more passive one. The comfortable seating system allows the user to alter the seat position, ensuring maximum comfort at all times. This combined with good driving characteristics makes **Azalea**® the ideal choice. **Azalea**® is a self propelling wheelchair with large rear wheels, and is complimented with an alternative transit version with small rear wheels called Rea™ **Azalea**® Assist.



Confidence in posture solutions



➤ **Maximised seating comfort and adjustments as key features**

With a combination of comfortable seating, easy adjustable seat angle and backrest recline and the possibility for self propelling or pushed by an attendant plus all the size adjustments makes the Rea™ **Azalea**® the natural choice when selecting a passive wheelchair.

➤ **Innumerable possibilities**

With three different frame sizes, a range of seat widths from 37 to 59 cm can be achieved. The wheelchair has individual width adjustments in the armrests and leg rest attachments as standard.

➤ **Tilt & recline makes life comfortable**

The Rea™ **Azalea**® and Rea™ **Azalea**® **Assist** have a 21° or 26° seat tilt respectively and a 30° backrest recline. Both tilt and recline can be easily adjusted manually or electrically by the assistant and also by the user for the electrical version. With the increased seat tilt for the Rea™ **Azalea**® **Assist** an enhanced body positioning and pressure relief can be achieved for maximum comfort and rest. The backrest has a biomechanical positioning hinge, which follows the user's natural back movement when adjusting the backrest angle.

➤ **Compatibility between product lines**

In order to keep the compatibility for the new Rea™ **Azalea**® with Rea™ **Bellis**, Rea™ **Assist** and Rea™ **Comfort** we use the existing Invacare® **Flex2** seating system as well as the same range of leg rest systems.





**Invacare® Flex2 assortment**  
Interchangeable **Flex2** seat system, with optional seat solutions.




**Length compensation**  
The leg rests and backrest include length compensation for preserving excellent support for the head and legs through height adjustable legs and backrests.




**Sliding balance point**  
When tilting the seat unit, the Rea™ **Azalea**® keeps its stability as well as its good driving performance.



**Minimized lifting at the front**  
The seat tilt mechanism is designed to give a minimum increased height at the front of the seat when tilted allowing the user to get under a table.

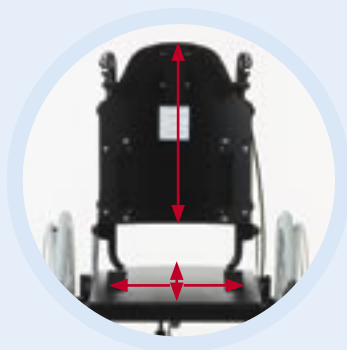

**Seat size**  
The size of the seat is adjustable in width, depth and backrest height to meet individual needs.

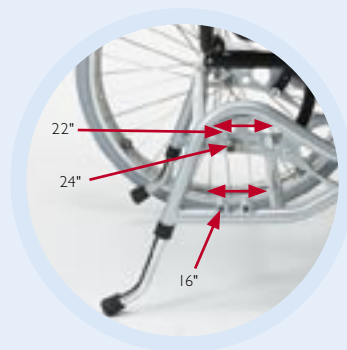

**Rear wheel adjustment**  
The seat height and the balance can be adjusted individually.




**Height adjustable padded arm- and side support**

A soft unit which supports the arms and the hips for increased comfort.



**Swing away table**


This type of table can easily be moved aside by a swing away mechanism and facilitating the transfer in and out in the chair.

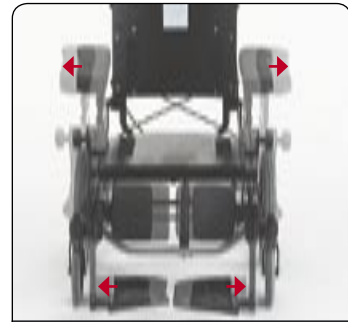

**Step less width adjustments**

Individual width adjustment in the armrests and leg rests.

**Technical data**

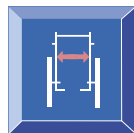


39-45/44-50/49-55 cm  
(39-49/44-54/49-59 cm)




43 - 50 cm



40/45 cm  
(45 cm)

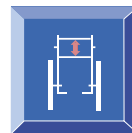

50 - 80 cm  
(without seat cushion)

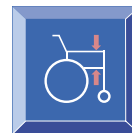

24-36 cm




33-50 cm



-1°-20°  
(-1°-25°)



0°-30°

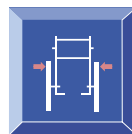


SW + 22 cm  
(SW + 23 cm)



90-125 cm



95-102 cm

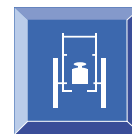

34 kg  
(SW 44 kg)



max. 135 kg



Transport weight  
16,5 kg



Crashtested  
ISO 7176-19\*

\* Our wheelchairs comply above ISO standards and with ISO 10542, and have additional features providing increased levels of occupant security and safety whilst travelling in a motor vehicle. Wheelchair users should however transfer to the vehicle seat and use the vehicle installed restraint system whenever it is feasible.

**Frame colour**



Pearl grey

**Upholstery**



Grey Dartex TR23



Grey plush TR 32

All information quoted is believed to be correct at time of print.  
Invacare® reserves the right to alter product specifications without prior consultation.  
Rea, Rea design and DSS (Dual Stability System) design are registered trademarks of Invacare® International.



INVACARE® Rea AB  
Box 200 · Växjövägen 303  
SE-343 75 DIÖ · Sweden  
Tel: +46 476 535 00 · Fax + 46 476 535 99  
www.invacare-rea.com

