



Yes, you can.[®]

rea Spirea^{®3}

An easily transportable wheelchair for everyday use

The Rea™ **Spirea^{®3}** is a durable, folding-type wheelchair that is simple to store away. It offers basic adjustments necessary for the user to sit steadily, whilst remaining a wheelchair that is easy to manoeuvre. **Spirea^{®3}** is easy to use, yet meets all the requirements of the user; assistant, therapist or technician. The **Spirea^{®3}** has been developed with our customers' needs in mind. Due to its high technical quality, the **Spirea^{®3}** is simple to disassemble, has low maintenance requirements and only a few simple, essential options.



Confidence in posture solutions



➤ **Simple to steer**

Rea™ Spirea^{®3} is designed to be steered by an assistant or by the user. The steering handles located at the back are adjustable in height and are available as accessories. The driving wheels are situated in an ideal position so that the wheelchair can be easily manoeuvred. The optimal steering characteristic of *Spirea^{®3}* is due to its stable framework, resulting in low rolling resistance.

➤ **Simple to manoeuvre**

Spirea^{®3} can be easily manoeuvred in narrow places. The striking armrests with handles for moving and handling can be easily removed. Similarly, supports can swing under the seat, for greater ease of access.

➤ **Simple to dissemble**

Spirea^{®3} is manufactured from high-quality durable aluminium, which provides excellent stability and low maintenance requirements. Most *Spirea^{®3}* parts are the same on the right side as on the left. The advantage of this is storage of fewer parts, which results in fewer part numbers in inventory.

Spirea^{®3} can be dissembled with three tools and assembled with just one. *Spirea^{®3}* has a unique chassis and various parts and components, including the wheels, armrests and leg supports can be used on other *Invacare® Rea™* wheelchairs.

➤ **Simple to transport**

Spirea^{®3} with a weight of 17,3 kg in standard design can be easily transported. It is simple to fold and carry.



Features and accessories



Rea™ Spirea^{®3} low seat height
With 22" flywheels Spirea^{®3} has a low seat height, and becomes a good kick-type wheelchair with free space below.



Adjustable backrest
Spirea^{®3} has a backrest which can be adjusted in angle and height. This makes it possible to form the backrest according to any individual needs.




Plaster cast leg support
Spirea^{®3} with plaster cast leg support and a fixed 80° leg support makes it extra stable.



Leg supports
Leg supports with fixed or adjustable footrest, fixed angle-adjustable plaster and amputation leg supports.



Armrests
Armrest with short, long or broad padding. Replaceable, detachable, hemiplegia armrests.

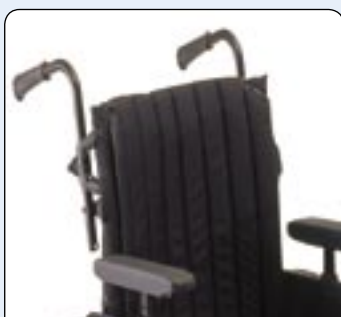


Steering handle located at the back
Height can be adjusted as desired by loosening the lever.



Safety against tipping
Height-adjustable. Combined with safety against toppling and trampling.



Stable Rea™ Spirea^{®3}
If a more stable wheelchair is required, attachment to the driving wheel located at the back can be used. Both alternative attachments to the driving wheel can be adjusted in height in order to obtain any sitting height.

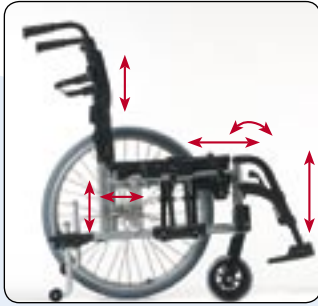

What can be adjusted

Back height, seat depth, angle of the flywheel, back angle, balancing of the driving wheels, as well as camber.




Standard design

Depth-adjustable feature, the backrest is adjustable in height and angle and it has adjustable limb feature and fixed steering handles. Leg supports and armrests are adjustable in height, as well as detachable. 24" compact driving wheel has 3° camber and has alternative balance and height adjustment. 6" compact flywheels, brake as well as safety against toppling are also incorporated.

Accessories

- Steering hoop
- Shape-adjustable seat with velcro strip
- Angle-adjustable leg supports
- Plaster case or amputation leg supports
- Angle- and depth-adjustable footrest
- Heel band
- Replaceable armrests
- Hemiplegia armrests
- Short or long padding strip of the armrest
- Broad and long padding of the armrest
- 120-200 mm flywheels; air/compact broad and narrow
- Drum brake
- Extended brake lever
- Unequal drive rings
- Spoke guards
- Reflex
- Table
- Crutch-holder
- Safety belt

Technical data

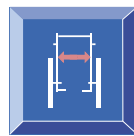


38, 40,5, 43, 45,5, 50,5 cm




40 cm (-5)
45 cm (-5)



38 - 50,5 cm

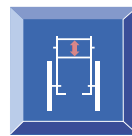

38 - 46 cm

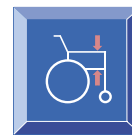

19 - 29 cm




32,7 - 52 cm




0° - 6°



± 12°

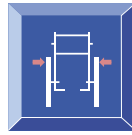


SB + 20 cm

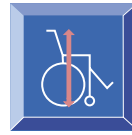


76 - 96,5 cm



90 - 115 cm

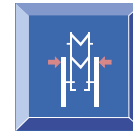

29 cm

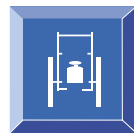

15,7 kg



Max 125 kg



Transportation
9,4 kg



Crash-tested*

* Invacares® wheel chairs fulfil the requirements of ISO 7176-19 and ISO 10542 which gives an increased degree of safety to the user and security when one travels by car. Therefore the wheelchair user can move up to the car seat and put the safety belt in the car.

Frame colour



Pearl grey

Upholstery



Black Jemima



Grey Plush

All information quoted is believed to be correct at time of print. Invacare® reserves the right to alter product specifications without prior consultation.



INVACARE® Rea AB
Box 200 · Växjövägen 303
SE-343 75 DIÖ · Sweden
Tel: +46 476 535 00 · Fax + 46 476 535 99
www.invacare-rea.com

